

Real-life performance of a novel antigen detection test on nasopharyngeal specimens for SARS-CoV-2 infection diagnosis: a prospective study

Authors:

L.Courtellemont¹, J. Guinard¹, C. Guillaume¹, S. Giaché², V. Rzepecki², A. Seve², C. Gubavu², K. Baud², H. Doize², C. Le Helloco², C. Lebegue², A. Tay², G.N. Cassuto³, G. Pialoux⁴, L. Hocqueloux², T. Prazuck²

1 Department of Virology, CHR Orléans, Orléans, France; 2 Department of Infectious and Tropical Diseases CHR Orleans, Orléans, France; 3 Laboratoires Drouot, Paris, France, 4 Department of Infectious and Tropical Diseases, CHU Tenon, Paris, France

Abstract

Introduction

The SARS-CoV-2 pandemic has become a major public health issue worldwide. Developing and evaluating rapid and easy-to-perform diagnostic tests is an absolute priority. The current prospective study was designed to assess diagnostic performances of an antigen-based rapid detection test (COVID-VIRO®) in a real-life setting.

Methods

Two nasopharyngeal specimens of symptomatic or asymptomatic adult patients hospitalized in the Infectious Diseases Department or voluntarily accessing the COVID-19 Screening Department of the Regional Hospital of Orléans, France, were concurrently collected. COVID VIRO® diagnostic specificity and sensitivity were assessed in comparison to real-time reverse transcriptase quantitative polymerase chain reaction (RT-qPCR) results. A subgroup of patients underwent an additional oropharyngeal and/or a saliva swab for rapid testing.

Results

121 patients already having a confirmed infection and 127 patients having no evidence of recent or ongoing infection were enrolled, for a total of 248 couple of nasopharyngeal swab specimens. Overall COVID-VIRO® sensitivity was 96.7% (IC: 93.5%-99.9%). In asymptomatic patients, patients having symptoms for more than 4 days and those having a RT-qPCR Cycle threshold value >32, sensitivity was of 100%, 95.8% and 96.9% respectively. The concordance

between RT-qPCR and COVID VIRO® rapid test was 100% for the 127 patients with no SARS-CoV-2 infection.

Conclusion

COVID-VIRO® test had 100% specificity and above 95% sensitivity, better than WHO recommendations (specificity ≥ 97 -100%, sensitivity ≥ 80 %). These rapid tests are particularly interesting for large-scale screening in Emergency Department, low resource settings and airports.

Introduction

At the end of 2019, a pneumonia of initially unknown origin was primarily reported to the World Health Organization (WHO) Country Office in China. On January 9th, 2020, the Chinese health authorities and the WHO announced the discovery of a novel coronavirus, firstly named 2019-nCoV, then officially SARS-CoV-2 (Severe Acute Respiratory Syndrome CoronaVirus 2). This virus, belonging to the coronavirus family but differing from SARS-CoV-1 and MERSCoV, is responsible for upper/lower respiratory tract infections known as COVID-19 (COronaVirus Disease 2019). COVID-19 incubation period is approximately 5.2 days and the most common onset symptoms are fever, cough, and fatigue (1). Since SARS-CoV-2 appeared in China, it has become a major public health issue worldwide. To date, more than 40 million cases have been detected worldwide (2) and the pandemic continues to spread unabated.

Minimizing testing delay seems to have the largest impact on reducing onward transmissions (3), and the availability of highly sensitive and specific tests is essential to quickly identify new cases and contain virus transmission.

Currently, the real-time reverse transcriptase quantitative polymerase chain reaction (RT-qPCR) assay is the gold-standard method to detect SARS-CoV-2 RNA in respiratory specimens such as nasopharyngeal or oropharyngeal swabs or broncho-alveolar lavage (4). However, performing RT-qPCR is expensive, time consuming, and requires special equipment and qualified operators. Faster, cheaper, and easier to use alternative tools could be represented by novel antigen-based rapid detection tests, or point of care tests (POCT) (5).

Recently, the WHO approved the first rapid detection test for a large-scale use in low- and middle-income countries (6), and French health regulation authorities authorized their utilization in medical settings (7).

Several different POCT have already been developed, with generally high specificity but variable sensitivity (8-13).

COVID-VIRO® (AAZ, Boulogne Billancourt, France) is one of the novel immunochromatographic tests designed to detect SARS-CoV-2 antigen in nasopharyngeal secretions within 15 min. Evaluating the diagnostic performance of COVID-VIRO® in the real life in comparison to the RT-qPCR as reference test is the principal aim of the current prospective study.

Materials and methods

Ethical approval

This study was approved by the Regional North West Ethics and Research Committee. A written informed consent was obtained from each participant.

Study population

People voluntarily accessing the COVID-19 Screening Department and SARS-CoV-2 positive patients hospitalized in the Infectious Diseases Department of the Centre Hospitalier Régional (CHR) of Orléans, France, or the Department of Infectious and Tropical Diseases of the Centre Hospitalier Universitaire (CHU) Tenon, Paris, France, from October 12th, 2020, to October 19th, 2020, were included in the study. Patient age was collected at inclusion, as well as symptom onset date for symptomatic patients. Suggestive symptoms were headache, fatigue, fever, or respiratory signs.

Participating patients underwent two concurrent nasopharyngeal swabs for RT-qPCR and COVID-VIRO® analysis, respectively.

Case definition

SARS-Cov2 positive subjects were either patients with a positive RT-qPCR at the time of the study sampling if done in parallel with the rapid test, either patients with a previous positive RT-qPCR within 5 days but a negative RT-qPCR at the time of study sampling.

SARS Cov2 negative subjects were patients with a negative RT-qPCR at the time of study sampling without any previous RT-qPCR.

Specimen collection

Two nasopharyngeal swabs were performed by trained personnel (nurse, doctors, or biologists). The collection of the two simultaneous samples was always carried out by the same operator.

A polyester-tipped flexible aluminum-shafted applicator (Microtest M4RT, Remel) was inserted into two of the nostrils until resistance was felt at the nasopharynx, then rotated 6 times and withdrawn. After swabbing, the swab applicator was cut off, and each absorbent swab was placed into a vial containing 3 mL of inactivating viral transport media. Half of nasopharyngeal swabs were immediately transferred to the Virology Unit of the CHR of Orléans hospital, Orléans, or Drouot Laboratory, Paris, to perform RT-qPCR, while the rapid antigen test was immediately on-site performed.

An additional oropharyngeal and/or saliva swab specimens were simultaneously collected in a subgroup of positive patients in order to determine the diagnostic reliability of these fluids in comparison to nasopharyngeal swab specimens. For oropharyngeal specimen, the swab was used to collect fluid on both sides of tonsillar arches and posterior pharynx. For saliva specimen, the sterile swab was used on the upper and lower gums all the way from the back to the front twice and immediately placed in the buffer vial.

Real-time RT-qPCR assays for the detection of SARS-CoV-2 RNA

Nucleic acid extraction was performed with automated Sample Preparation System MGISP-960 (MGI, China). Specific real-time RT-qPCR assays target three SARS-CoV-2 genes, namely ORF1ab, S and N genes (TaqPath Covid-19 Multiplex RT-PCR, Thermofisher). Genome amplification was performed using QuantStudio5 (Applied Biosystems). Results interpretation was performed according to manufacturer instructions. The assay includes an RNA internal extraction control and amplification control. Samples showing an exponential growth curve and a Cycle threshold (Ct) value < 37 were considered as positive. A Ct unique value > 37 was considered as negative.

Rapid antigen test

COVID-VIRO® (AAZ, Boulogne Billancourt, France) is a membrane-based immunochromatography assay detecting SARS-CoV-2 nucleocapsid antigen (N-protein) in nasopharyngeal samples through monoclonal antibodies. A second

monoclonal antibody is conjugate to colloidal gold particles which are captured in reaction membrane. The test was performed according to manufacturer instruction by mixing nasopharyngeal secretions with 300 μ L of dilution buffer in a tube. One minute, 4 drops were added in the appropriate well. When nasopharyngeal secretions cross the strip, a passive diffusion allows the solubilized conjugate to migrate with the sample and react with the anti-SARS-CoV-2 antibodies immobilized on the membrane. A control line allows assessing the correct migration of sample and the reliability of the test. Visual interpretation of results is performed 15 min after (Figure 1).

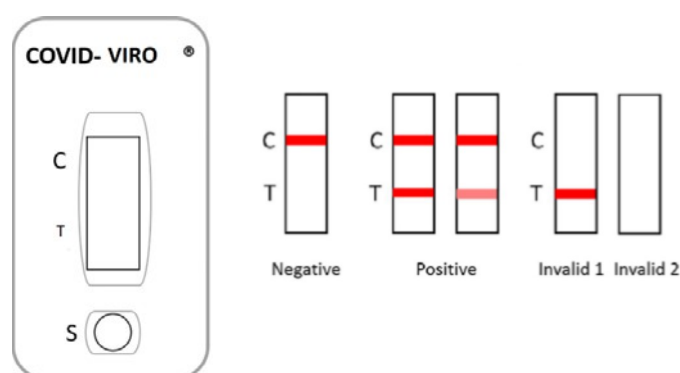


Figure 1: Interpretation of results for COVID-VIRO[®] (AAZ, LMB)

Data analysis

Population characteristics are reported as percentage, mean and median values, standard deviation, and range. Data was analyzed in the Infectious Diseases Department.

To determine the diagnostic value of COVID-VIRO[®], the study population was divided into two groups:

1. Already confirmed RT-qPCR positive patients. Comparison between RT-qPCR and COVID-VIRO[®] results in these patients was used to assess diagnostic test sensitivity.
2. Non-selected symptomatic or asymptomatic patients voluntarily accessing the COVID-19 Screening Department in order to detect a possible SARS-CoV-2 infection. Analysing data, RT-qPCR positive patients were added to the first group in order to assess diagnostic test sensitivity. Conversely, RT-qPCR negative patients were selected to measure the specificity of the rapid test.

COVID-VIRO® specificity and sensitivity were calculated using the RT-qPCR results as reference test, according to the following formulas:

$$\text{Specificity (\%)} = 100 \times [\text{Negative} / (\text{Negative} + \text{Positive})]$$

$$\text{Sensitivity (\%)} = 100 \times [\text{Positive} / (\text{Positive} + \text{Negative})]$$

The size of the study population was calculated on the basis of a 95% sensitivity with a lower margin of the confidence interval over 91% and a 99% specificity with a lower margin over 95% according to the WHO recommendations for antigenic rapid test. T Student test was used to compare means.

Results

121 patients with a SARS-CoV-2 confirmed infection and 127 patients with no evidence of SARS-CoV-2 recent or ongoing infection were enrolled in the study. A total of 248 couples of nasopharyngeal swabs were analyzed. Of these, 228 were collected in Orléans, and 20 in Paris.

The sex ratio of the study population was 0,9 (117 men and 131 women). The median and mean age was 38 and 43 years old, respectively (range: 18-96). Figure 2 shows the flow chart of included subjects in the study. One patient exhibiting only one positive target (gene S, Ct 36) at RT-qPCR analysis, was considered as negative according to the French Society of Microbiology criteria and the extraction kit instructions, and was included in the group of SARS-CoV2 negative patients.

Among the 121 SARS-CoV-2 confirmed infected patients, 17 patients were hospitalized (14%) and 99 patients (81.8%) were symptomatic. The median time of symptoms duration before sampling was 5 days (mean: 5.3 days, range: 1-20).

N gene mean Ct value was 29 (range: 15-34) in asymptomatic SARS-CoV-2 confirmed patients, compared to 27 (range: 13-35) in symptomatic positive patients. The difference was not significant ($p < 0.05$). In patients tested within 0 to 4 days from symptoms onset, the N gene mean Ct value was 25 (range: 13-35), versus 28 (range: 18-33) in those tested within 5 and 7 days, and 30 (range: 21-35) in those having symptoms for more than 7 days, showing a constant decrease in viral carriage. Although some false negative patients had higher Ct

value (i.e: low RNA carriage), COVID-VIRO® was able to detect the antigen in 15 patients having symptoms for more than 7 days.

Among the 121 RT-qPCR positive patients, 4 had a negative COVID-VIRO® result (3.3% false negative). The overall sensitivity of the POCT is estimated to be 96.7% (IC: 93.5%-99.9%) (table 1). There was no false negative result obtained with COVID-VIRO® among the 22 asymptomatic patients.

Table 2 shows the COVID-VIRO® performances according to the Ct value and the delay of symptoms onset. COVID-VIRO® sensitivity was extremely high among patients having a Gene N, S or ORF Ct-values > 32, considering that 32 out of 33 patients were tested positive (sensitivity: 96.9% (95% IC: 91.1% - 100%).

Table 3 shows the characteristics of the four COVID-VIRO® false negative cases. Three out of four had Ct values \geq 32 and were clinically considered as non-contagious.

Twenty positive patients having a previous positive RT-qPCR were tested negative when the second RT-qPCR was dually performed with the POCT. All these patients had a positive COVID-VIRO® result (mainly weak or very weak line) suggesting that it could be still positive some days after PCR negativization.

Among the 127 patients with no SARS-CoV2 infection, no false positive result was observed and the concordance between RT-qPCR and COVID-VIRO® was 100%. Therefore, COVID-VIRO® specificity is estimated to be 100%.

Additionally, 48 patients having a positive COVID-VIRO® test on nasopharyngeal swab specimen accepted undergoing a simultaneous oropharyngeal (34 patients) or saliva swab (14 patients). COVID-VIRO® was positive in 24 out of 34 and 0 out of 14 patients on oropharyngeal and saliva specimens, respectively. Sensitivity was 70.6% for oropharyngeal specimens and 0% for saliva specimens.

Discussion

This prospective observational study aims to evaluate the performance of a POCT designed to detect SARS-CoV-2 antigen from a nasopharyngeal swab directly after sampling and giving the result within 15 minutes. The diagnostic value of COVID-VIRO® was determined using RT-qPCR as gold standard in a real-life community setting. In our study, the test sensitivity reached 96.7% and specificity was found to be 100%, with no false positive observed. Although if visual interpretation of the result is performed after 15 min, almost always a positive result appeared within the first five minutes, not rarely within one minute.

To date, few studies have already evaluated the performance of such POCT to detect SARS-CoV-2 antigen on nasopharyngeal fresh swab in the real life. The only prospective studies available concern the Panbio COVID-19 Ag Rapid test (Abbott). While the manufacturer reported a high sensitivity (93.3%; 95 CI 83.8-98.2) in a high endemic setting in Brazil (14), other independent cohort studies did not show such performances. In 257 symptomatic and asymptomatic patients enrolled at the Emergency Department and Primary Health Care Setting in Spain, overall sensitivity was 73.3%, reaching 86.5% among patients having symptoms for less than seven days (15). In another multicentric study performed on 200 COVID-19 RT-PCR positive patients, POCT sensitivity was 72.6% (95% CI: 64.5-79.9%) in the Netherlands and 81.0% (95% CI: 69.0–89.8%) in Aruba. Test sensitivity was as high as 95.2% (95% CI: 89.3-98.5%) in patient with RT-qPCR test positivity to Ct-values < 32. (16).

Even in our study, patients were stratified by the Ct value to evaluate the test sensitivity at different level of nasopharyngeal viral load. COVID-VIRO® sensitivity remained extremely high even when Ct values were >32 (96.9%; 95% IC: 91.1% -100%).

Detection of viral RNA on nasopharyngeal samples is not necessarily linked to infectiousness (17). Several factors determine viral transmission risk: these include whether a virus is still viable, the amount of replicative virus estimated by the Ct, the presence of respiratory symptoms, the individual's local mucosal immune response to the virus, and the behaviour of the infected individual and their contacts (18). However, in the present study the number of viral particles estimated by the Ct value did not differ in asymptomatic and symptomatic SARS-CoV-2 infected patients. COVID-VIRO® appears to be as sensitive as RT-qPCR to detect infected patients in a limited number of asymptomatic patients.

Available data report that RNA viral load rapidly decreases after the onset of symptoms, and infectiousness generally declines within 7–10 days (17, 19-22). Stratifying our patients by symptoms duration at sampling, we observed a similar decline from day 1 to day 14 in the mean Ct value enregistered. Just one hospitalized patient aged of 90's, exhibited a Ct value of S :21, N: 22, ORF: 21 at day 10. Anyway, our POCT was able to detect the antigen as well as within than after 4 days of symptom onset.

Considering discordant results, the analytical performances depend on different factors including the viral load, the quality of the specimen and the modalities of processing. The two nasopharyngeal swabs were concurrently performed by the same operator, but we can hypothesize that a greater quantity of secretions and therefore of virus is concentrated on the first swab instead of on the second. Unfortunately, we do not know the temporal order in which each swab was performed, due to the methodology used, and waiting one hour between the two samples was infeasible.

Analysing POCT results on oropharyngeal and saliva specimens of positive patients, we quickly realized that sensitivity dramatically drops. Even if nasopharyngeal swab is not comfortable for patients, this specimen should be privileged for the diagnosis of SARS-CoV-2 infection, whether obtained through PCR or POCT.

This study has several limitations. First, the date of symptoms onset was reported by patients and may not always be accurate, leading to an inaccurate stratification of patients. Second, the number of asymptomatic patients is rather limited to obtain conclusive data, even if we did not observe any discordance between tests in this group of patients.

In France, nurses, pharmaceuticals, and general practitioners have recently been authorized to perform POCT in medical settings (7). Data obtained from the current study could reassure health authorities in expanding access to POCT that, in addition to being quick and easy to use, are also reliable. Currently, symptomatic people risk spending some days waiting for the execution and then for the result of a SARS-CoV-2 RT-qPCR screening test. During this time, unaware infected people waiting for the result are more at risk to transmit the virus than after having the official result of their test. A more rapid diagnosis and the subsequent contact tracing would certainly positively impact on the containment of transmission.

In the Emergency Departments, POCT could also be used to quickly recognise asymptomatic positive from negative patients, avoiding SARS-CoV-2 nosocomial infection. Furthermore, such tests would likely be useful in low- and middle-income countries and at airports to limit the worldwide spread of SARS-CoV-2.

Conclusion

Performance of COVID-VIRO® (AAZ, Boulogne Billancourt, France) was a quite a reliable test for SARS-CoV-2 diagnosis. The sensitivity and specificity of the evaluated POCT in nasopharyngeal swabs were 96.7 % (95% CI 9 – 99.9%) and 100% respectively. To date, this is the unique COVID-19 antigenic rapid test fulfilling the WHO's recommendations for a screening test (sensitivity $\geq 80\%$, specificity $\geq 97-100\%$). Unfortunately, performing the test on oropharyngeal or salivary samples is not enough reliable.

Bibliography

1. Rothan HA, Byrareddy SN. The epidemiology and pathogenesis of coronavirus disease (COVID-19) outbreak. *Journal of Autoimmunity*. 2020; 109: 102433. <https://doi.org/10.1016/j.jaut.2020.102433> PMID: 32113704
2. Dong E, Du H, Gardner L. An interactive web-based dashboard to track COVID-19 in real time. *The Lancet Infectious Diseases*. 2020
3. Kretzschmar ME, Rozhnova G, Bootsma MCJ, van Boven M, van de Wijkstra JHHM, Bonten MJM. Impact of delays on effectiveness of contact tracing strategies for COVID-19: a modelling study. *Lancet Public Heal*. 2020
4. Gorbalenya AE, Baker SC, Baric RS, de Groot RJ, Drosten C, Gulyaeva AA, et al. Severe acute respiratory syndrome-related coronavirus: The species and its viruses—a statement of the Coronavirus Study Group. *Microbiology*; 2020 Feb. <https://doi.org/10.1101/2020.02.07.937862>
5. Dinnes J, Deeks JJ, Adriano A, et al. Rapid, point-of-care antigen and molecular-based tests for diagnosis of SARS-CoV-2 infection. *Cochrane*

- Database Syst Rev. 2020;8:CD013705. Published 2020 Aug 26. doi:10.1002/14651858.CD013705
6. WHO.com. retrieved October 26 2020 [Internet]. Available from: <https://www.who.int/news/item/28-09-2020-global-partnership-to-make-available-120-million-affordable-quality-covid-19-rapid-tests-for-low--and-middle-income-countries>
 7. Arrêté du 16 octobre 2020 modifiant l'arrêté du 10 juillet 2020 prescrivant les mesures générales nécessaires pour faire face à l'épidémie de covid-19 dans les territoires sortis de l'état d'urgence sanitaire et dans ceux où il a été prorogé – Légifrance JORF n°0253 du 17 octobre 2020
 8. Finddx.org. retrieved October 13 2020 [Internet]. Available from: <https://www.finddx.org/covid-19/pipeline/>
 9. Hirotsu Y, Maejima M, Shibusawa M, Nagakubo Y, Hosaka K, Amemiya K, et al. Comparison of automated SARS-CoV-2 antigen test for COVID-19 infection with quantitative RT-PCR using 313 nasopharyngeal swabs, including from seven serially followed patients. *Int J Infect Dis*. 2020;
 10. Scohy A, Anantharajah A, Bodéus M, Kabamba-Mukadi B, Verroken A, RodriguezVillalobos H. Low performance of rapid antigen detection test as frontline testing for COVID-19 diagnosis. *J Clin Virol*. 2020;
 11. Guglielmi G. Fast coronavirus tests: what they can and can't do. *Nature*. 2020.
 12. Linares M, Pérez Tanoira R, Romanyk J, Pérez García F, Gómez-Herruz P, Arroyo T, et al. Panbio antigen rapid test is reliable to diagnose SARS-CoV-2 infection in the first 7 days after the onset of symptoms. *medRxiv*. 2020;
 13. Dinnes J, Deeks JJ, Adriano A, Berhane S, Davenport C, Dittrich S, et al. Rapid, point-of-care antigen and molecular-based tests for diagnosis of SARS-CoV-2 infection. *Cochrane Database Syst Rev*. 2020.
 14. globalpointofcare.abbott. retrieved October 13 2020.
 15. Linares M, Pérez Tanoira R, Romanyk J, Pérez García F, Gómez-Herruz P, Arroyo T, et al. Panbio antigen rapid test is reliable to diagnose SARS-CoV-2 infection in the first 7 days after the onset of symptoms. *medRxiv*. 2020;
 16. Gremmels H, Inkel BMF, Schuurman R., Ritger NAM Rodriguez O. et al. Real-life validation of the Panbio COVID-19 Antigen Rapid Test (Abbott) in community dwelling subjects with symptoms of potential SARS-CoV-2 infection. *medRxiv*. 2020;

17. Walsh KA , Jordan K , Clyne B , Rohde D , Drummond L , Byrne P , et al. SARS-CoV-2 detection, viral load and infectivity over the course of an infection. *J Infection* 2020; **81** (3):357–71.
18. World Health Organization. Criteria for releasing COVID-19 patients from isolation 2020 [updated 17 Jun 2020; cited 2020 27 Aug]. Available from: <https://www.who.int/publications/i/item/criteria-for-releasing-covid-19-patients-from-isolation>.
19. Cevik M , Tate M , Lloyd O , Maraolo AE , Schafers J , Ho A . *SARS-CoV-2, SARS-CoV-1 and MERS-CoV viral load dynamics, duration of viral shedding and infectiousness: a living systematic review and meta-analysis*, 2020. medRxiv; 2020. 20162107, 07.25.
20. Jefferson T, Spencer E, Brassey J, Heneghan C. *Viral cultures for COVID-19 infectivity assessment. Systematic review*, 2020. medRxiv; 2020. 20167932, 08.04.
21. Park M, Pawliuk C , Nguyen T , Griffitt A , Dix-Cooper L , Fourik N , et al. *Determining the period of communicability of SARS-CoV-2: a rapid review of the literature*, 2020. medRxiv; 2020. 20163873, 07.28.
22. Walsh K.A., Spillane S., Comber L., Cardwell K. et al. The duration of infectiousness of individuals infected with SARS-CoV-2. *J Infection* 2020, <https://doi.org/10.1016/j.jinf.2020.10.009>

| | | COVID-VIRO® | | COVID-VIRO® performances | |
|----------------------------------|-----------------|-------------|----------|---------------------------|-----------------|
| | | Positive | Negative | Sensitivity(%, 95% IC) | Specificity (%) |
| All patients | Positive RT-PCR | 113 | 4 | 96.7% (93.3-99.9%) | |
| | Negative RT-PCR | 0 | 109 | | 100% |
| Asymptomatic COVID-19 + patients | Positive RT-PCR | 22 | 0 | 100% | |

Table 1: Performances of the COVID-VIRO® antigenic rapid test in the overall population and in the group of asymptomatic SARS-CoV2 infected patients

| | | COVID-VIRO® | | Sensitivity (% , 95% IC) |
|--|-----------------|-------------|----------|--------------------------|
| | | Positive | Negative | |
| Patients with high viral carriage (Ct≤28) | Positive RT-PCR | 66 | 1 | 98.5% (95.6-100%) |
| Patients with high and moderate viral carriage (Ct<32) | Positive RT-PCR | 88 | 3 | 96.7% (93.0-100%) |
| Delay from symptoms onset≤4 days | Positive RT-PCR | 36* | 2 | 94,7 % (87.6-100%) |
| Delay from symptom onset >4 days | Positive RT-PCR | 62* | 2 | 95.8 % (90.2-100%) |

Table 2: Sensitivity of the COVID-VIRO® antigenic rapid test in comparison to RT-PCR according to viral carriage and delay from symptom onset

| Patient | Age range | Sex | Symptoms | S, N, ORF Gene Ct respectively | Delay from symptom onset |
|---------|-----------|-----|----------|--------------------------------|--------------------------|
| #1 | 20-30 | F | Yes | 31, 34, 31 | 5 |
| #2 | 20-30 | M | Yes | 32, 32, 31 | 1 |
| #3 | 60-70 | F | Yes | Neg, 32, 32 | 16 |
| #4 | 40-50 | M | Yes | 32, 28, 28 | 3 |

Table 3: Characteristics of the four discordant positives RT-PCR negative COVID-VIRO®. None of patient was hospitalized.

Figure 2: Flow Chart (*: one patient re-assigned, RT-PCR Ct>36 one target alone, false positive RT-qPCR)

